

*Sleek, modern and stylish, Phoenix Orthodontics premier line of MIM (metal injection molded) buccal tubes incorporate the latest design technology with a growing industry demand for comfortable, easy-to-use appliances that simplify treatment and reduce valuable chair time.*

**ORION LP**  
low profile buccal tubes

\$3.95 ea  
bondable




**Unique Direct Bond Pad Design**

ORIONLP features a compound radius on the pad that mirrors the unique contour of a molar's clinical crown.

A uniquely crafted distal gingival cutout enables bonding of partially erupted molars without sacrificing the size and strength of the pad.

Choose from extended indent pad or small indent pad.

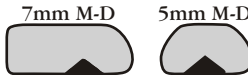

**80 Gauge Foil Mesh**

ORIONLP bondable buccal tubes feature 80 gauge foil mesh direct bond pads for increased bond strength.



Convertible upper & lower  
1st molar. Roth and MBT Rx.

\$4.50 ea bondable

**Maximum Rotational Control**

At a full 4.2mm, ORIONLP Tubes allow for greater rotational control than smaller buccal tubes.

**Distal Notch**

In order to meet the need for ample anchorage, ORIONLP features a distal notch for ligation and tie-back of auxiliaries.

**Fully Funneled Entrance**

One of the most appealing features of the ORIONLP is its trumpeted entrance to facilitate easy wire insertion.



ORIONLP on welding flange

Specify WELDABLE or BONDABLE

**MAXILLARY 1st & 2nd Molar**    **MANDIBULAR 1st & 2nd Molar**

	<u>torque</u>	<u>offset</u>		<u>torque</u>	<u>offset</u>
Roth	-10	8	Roth	-25	-5
MBT	-14	10	MBT	-20	0
Std.	0	0	Std.	0	0

- 80 gauge foil mesh for superior bond strength.
- Compound-contour pad w/ buccal indent for snug fit.
- Pad specs: 5mm M-D / 3mm height.
- Slim, malleable hook allows easy engagement.
- Smooth, rounded tube contours = patient comfort.
- Distal gingival cutout lessens tissue contact.

**FE-B7**

(fully erupted 2nd molar bondable tube)



\$3.95 ea  
bondable

MAXILLARY (MBT / Roth)

<u>torque</u>	<u>tip</u>	<u>offset</u>
-14	0	10

MANDIBULAR (MBT)

<u>torque</u>	<u>tip</u>	<u>offset</u>
-14	0	0

MANDIBULAR (Roth)

<u>torque</u>	<u>tip</u>	<u>offset</u>
-30	-1	4

**EASY PLACEMENT.** Facial surface indent on buccal tube assists in bonding and positioning for precise mesio-buccal fit.

**EASY FIT.** Elongated, curved, mesially extended contoured base (3.5mm M-D) reduces unwanted bond failure.

**EASY WIRE INSERTION.** Funneled entrance effectively guides wire into the buccal tube.

**EASY ON YOUR PATIENT.** Smooth, rounded low profile buccal tube (2mm M-D) provides maximum patient comfort. Comfortable hook design tilts away from the gingiva to avoid tissue impingement.



\$3.50 ea  
bondable

MAXILLARY torque	MANDIBULAR torque
0	0
-10	-12
-10	-12
	-25
	-25

