

PUMICE, ADHESIVE REMOVAL & FINISHING SYSTEM

**1st & Final
 ORTHODONTIC PUMICE**



PRE-PACKAGED:

200 sterile cups for individual patient use. Pre-mixed-no messy mixing with water. Just open and apply with prophyl cup.

IMPROVED BOND STRENGTH:

University tests show that 1st & Final will pumice more effectively and rinse cleaner for increased bond strength.

NO ADDED CHEMICALS:

No flavoring, fluoride or oils which could leave a film that may inhibit or weaken the bond.

ENAMEL POLISHER:

Can be used as the final polishing step with the *Renew Finishing System*.

RENEW FINISHING SYSTEM BURS:

The *118S* is a special fluted bur designed for controlled removal of bulk adhesive after the brackets have been removed. The *118L* is a longer shank version of the *118S* with the same excellent adhesive removal properties. These non-tapered friction grip burs cut quickly and do not provoke fissures. The *218* bur is an 18 fluted design that is tapered to reach into those tight areas along the gingiva and occlusal when removing residual adhesive. The wide flutes self clean. These burs are the first step in the *Renew Finishing System*.



RENEW SYSTEM BUR 118S (RSB)




RENEW SYSTEM BUR 118L (RSBL)



RENEW SYSTEM BUR 218 (RSB2)

RENEW SYSTEM POINTS:

These points will remove small residual amounts of adhesive from the tooth surface and will polish the enamel without damage. It's rubberized tapered shape facilitates removal of adhesive close to the gingiva. *383 Point* is the second step in the *Renew Finishing System*.

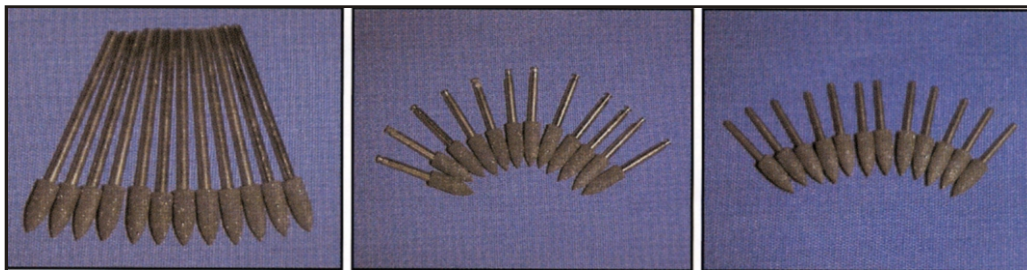

Straight Handpiece (RSPHP)

Latch (RSPL)

Friction Grip (RSP)

DIAMOND BUR 129:

A medium abrasion diamond bur recommended for roughening and preparing plastic, composite restorations and metal surfaces. An excellent quality medium bur at an affordable price.



DIAMOND BUR 129 (db)